

NEBULIZER




MOUTHPIECE

A nebulizer turns medicine into a mist that you can breathe deeply into your airways. Nebulizers are used when a person with asthma has trouble using a metered dose inhaler. Some medicines can be mixed together in a nebulizer - check with your doctor.




1: PREP MEDS

Put the medicine in the nebulizer cup and attach the cap.

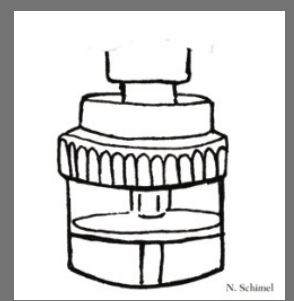


2: JOIN PIECES

Put T-piece on top of nebulizer cup. Attach the mouthpiece on one end of the T-piece and the wrinkled tube at the other end of the T-piece.

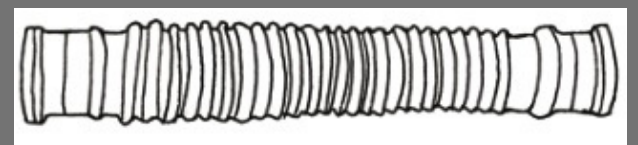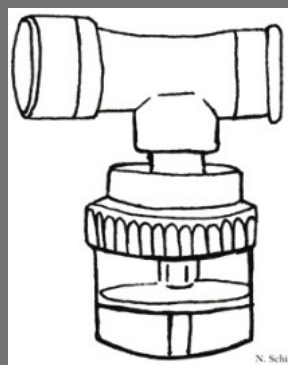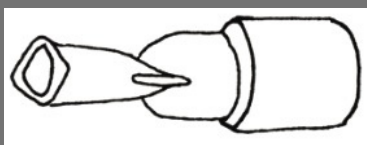

3: ATTACH TUBES

Attach one end of the clear, skinny tube to the bottom of the nebulizer cup and the other end to the nebulizer machine.


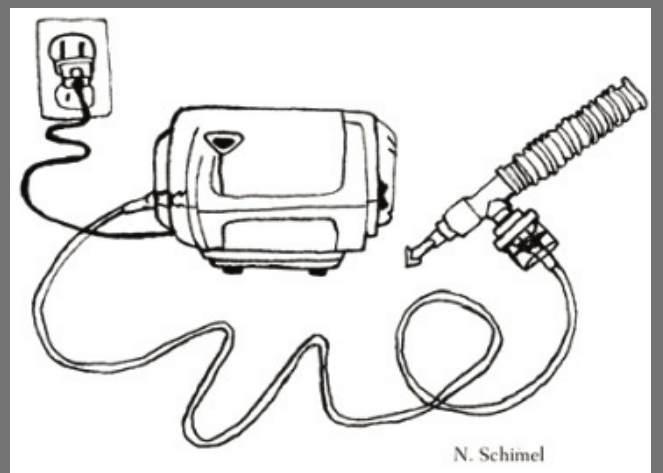


4: GET STARTED

Have the person sit straight up.

Put the mouthpiece between the teeth and seal the lips.

Turn the machine on.

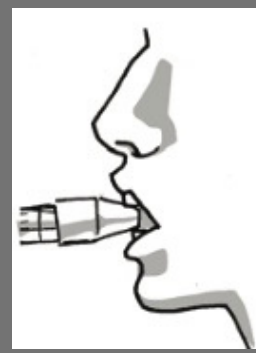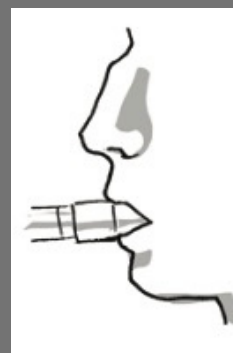

Breathe in slowly through the mouth, then exhale through the mouthpiece. The treatment is done when all the medicine is gone and there is no more mist. Tap the side of the medicine cup to force all the medicine through.

AFTER USING A NEBULIZER + MOUTHPIECE

RINSE AND SPIT



CLEAN & DISINFECT

Do not soak, wash or rinse the skinny tubing or the machine!

Clean after every use: Wash the medicine cup, t-piece, and wrinkled tube in soapy water.
Rinse with water and air dry

Disinfect monthly: Mix 1 cup vinegar and 1 cup water.
Soak the medicine cup, t-piece, and wrinkled tube overnight.
Rinse with water and air dry.

CARE FOR THE MACHINE

Wipe with a damp cloth as needed

Check the filter every month and change it when it becomes discolored.

Send the machine to be serviced every 5 years.